
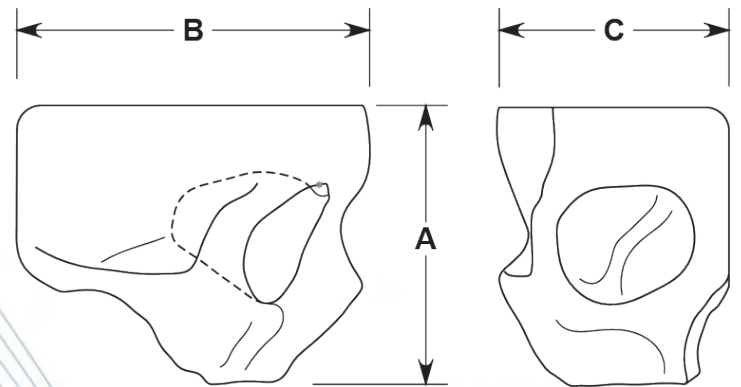


Surgical Implants

- Designed to restore non-load bearing bony structures of the orbit
- May be modified with a scalpel, scissors, or burr to fit a variety of defects
- Fully porous material for tissue ingrowth
- For major shape adjustments, submerge implant in a hot, sterile saline bath (82°-100°C, 180°-212°F) and allow to cool in desired shape
- Fixation screws or sutures can be placed in any desired location
- Individually packaged, MR Safe, Sterile, 5 Year Shelf Life

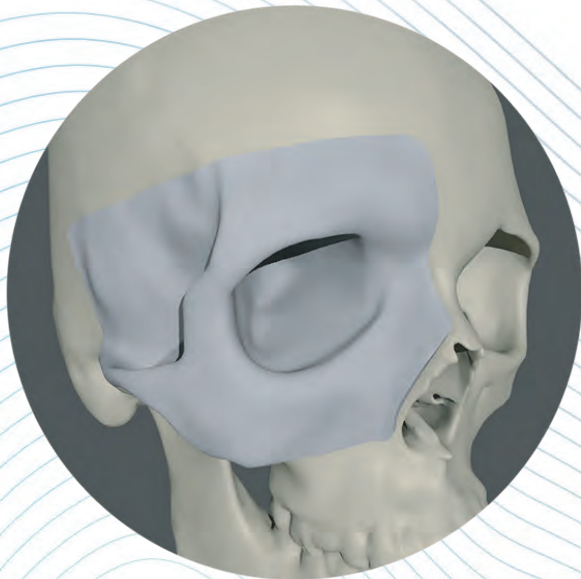

CAT#	Description	A	B	C
4226	Left	77mm	97mm	63mm
4227	Right	77mm	97mm	63mm