



- This implant design provides an excellent option for closure following the supraorbital approach.
- Incorporates an impervious layer on the dura facing side
- Scissors can be used to modify the shape of the implant approximately 5mm past the defect margins, providing a flange for direct fixation to the underlying bone.
- Typical fixation utilizes 5mm long self-drilling screws (not included).
- A scalpel or burr can be utilized to smooth the transition to the patient's bone directly in vivo.
- Retention Shape Memory™ provides malleability by hand.
- Provided sterile by ETO
- MR Safe, 5 Year Shelf Life

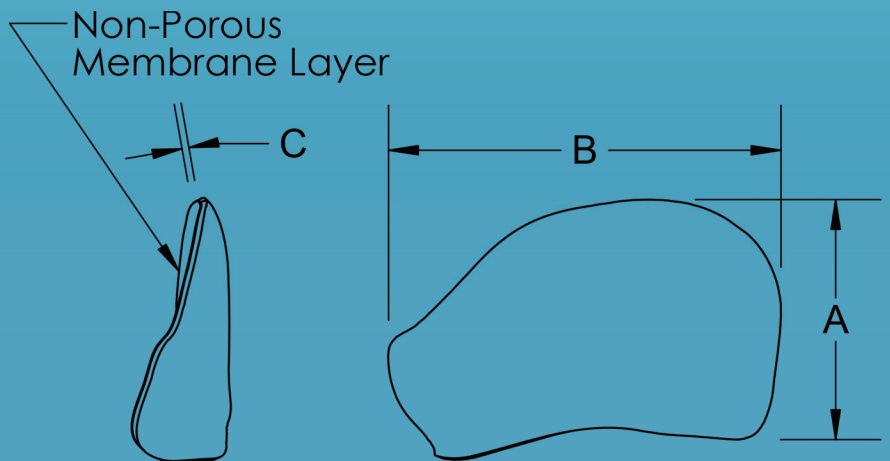

Cat #	Description	A	B	C
4572	Supraorbital Eyebrow Closure - Left	34mm	55mm	1.5mm
4573	Supraorbital Eyebrow Closure - Right	34mm	55mm	1.5mm

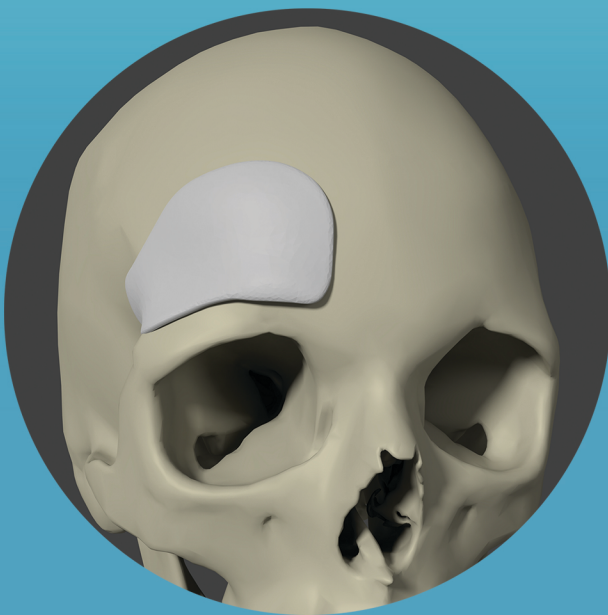